

# **Part III Patient Encounters and Performance Skills (PEPS®) Examination**

**SITE INFORMATION & EQUIPMENT LIST**

The information listed in this Part III PEPS Examination Site Information and Equipment List is for informational purposes only. There may be additional equipment or supplies on-site during the exam that may not be reflected in this manual. Additionally, equipment and supplies may change without warning.

This information is not intended to prepare candidates for the exam, but to provide a general overview of equipment and supplies used during the exam. Candidates are reminded to thoroughly review the evaluation forms and Candidate Guide to adequately prepare for the Part III PEPS Examination.

## **Part III PEPS EQUIPMENT QUICK REFERENCE:**

### **All Stations:**

- Laptop device (specific information will be available February 2024)

### **Patient Encounter Stations:**

- Reliance 6200-L Chair

### **Anterior Segment Station:**


- Reliance 6200-L Chair
- 7900-IC Floor Stand
- Haag-Streit 900 BI Slit Lamp Video System
- Goldmann-style applanation tonometer with tonometer prism
- Volk 4-Mirror Glass Gonio Fundus Lens with Flange

### **Posterior Segment Station:**

- Reliance 6200-L Chair
- 7900-IC Floor Stand
- Haag-Streit 900 BI Slit Lamp Video System
- “Modified” Keeler BIO Vantage Plus Slimline Wireless Plus Indirect System with +2.00 D Oculars
- Volk 20D Condensing Lens
- Volk VDGTL1 Digital High Mag Imaging Lens
- Volk Digital Wide Field Lens
- Volk 78D & 90D Lenses

## CANDIDATE DEVICE

Candidate must log into the device during check-in. Candidate will have the device during orientation and through the exam. Please review separate user guide for specific instructions on using the platform.

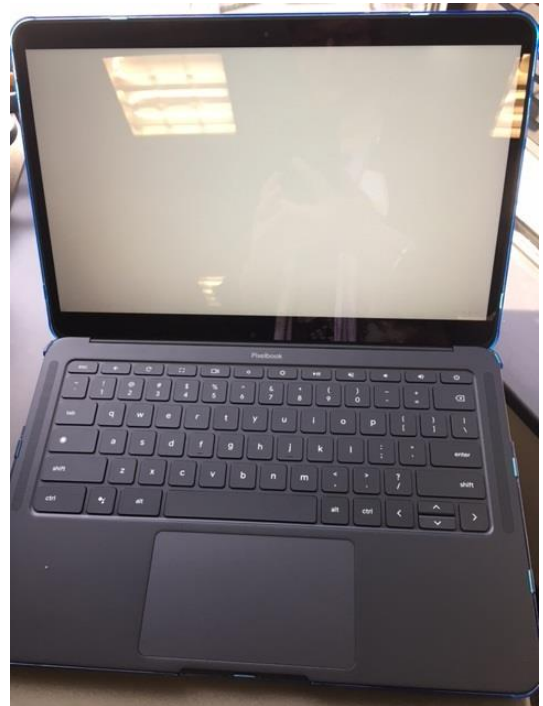

10 minutes before the exam is scheduled to start, the candidate's device will no longer be able to access orientation information.

During the exam, the Candidate will be able to use their device to view/submit:

- Skill station and encounter station overview information
- View patient encounter clinical data and tests
- Complete evaluation notes

## ALL STATIONS

**The following information is identical in all exam rooms**

### **MAIN OVERHEAD LIGHT SWITCH, DIMMER LIGHT SWITCH, AND SLIDER**


Located behind the computer workstation

Location of this switch is the same in EVERY exam room

Main overhead light switch is located on the LEFT

Dimmer light switch is located on the RIGHT

Slider to adjust dimmer brightness is on the RIGHT side of the dimmer switch

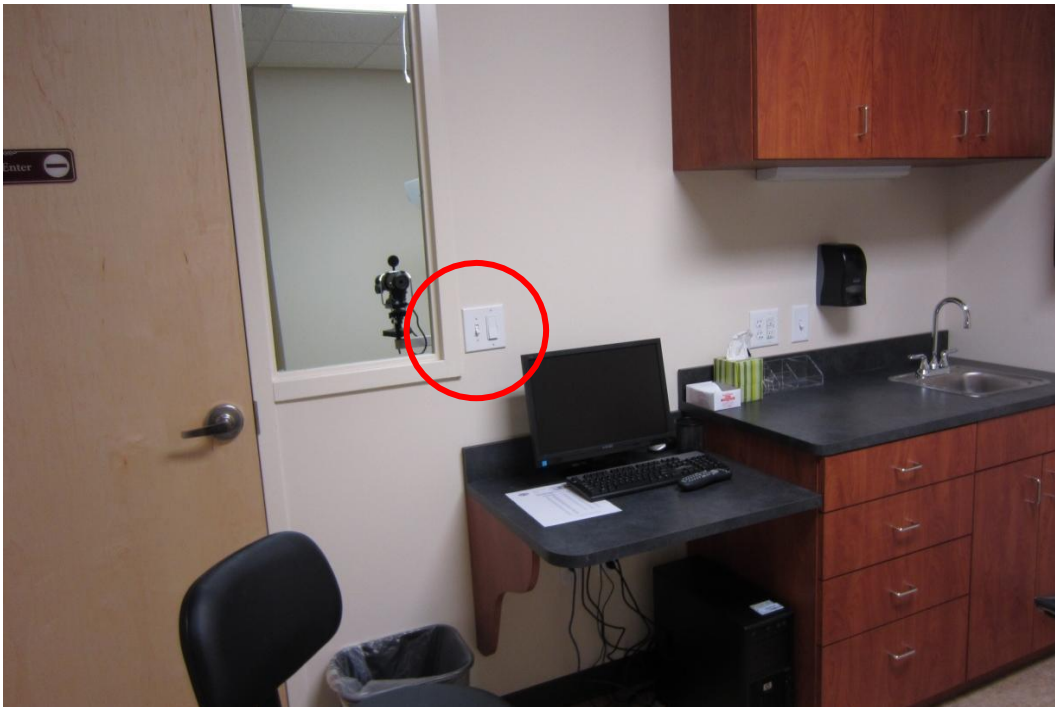

Lighting shown on this page for skills rooms

**Normal room illumination is when the main overhead lights are on**



**Slider to adjust dimmer brightness**

## ALL STATIONS

**The following information is identical in all exam rooms**

### UNDER-CABINET LIGHT SWITCH

Located above the sink counter, by soap dispenser

Location of this switch is the same in EVERY exam room



## ENCOUNTER STATIONS

### Patient Encounter Room Views





## ANTERIOR AND POSTERIOR EXAM CHAIR POSITIONING

### Exam Chair Overview Reliance 6200-L Chair



## Exam Chair Overview


### Reliance 6200-L Chair

#### To recline the exam chair:

Pull DOWN on the lever located on the back of the chair

Lower the chair and position appropriately

Release the lever while lowered to lock the chair in place

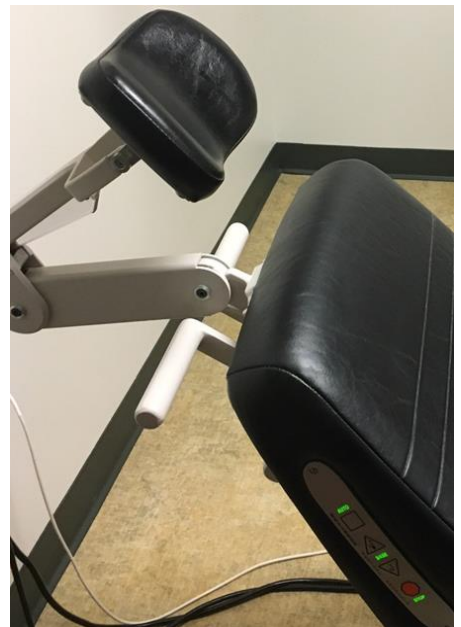

To rotate the chair, push the silver handle towards the floor. When the chair is in position, pull handle to the original position to lock in place.




## ANTERIOR AND POSTERIOR STATIONS SLIT LAMP OVERVIEW

### Slit Lamp


Haag-Streit 900 BI Slit  
Lamp Video System

**CAMERA IS ON THE LEFT OCULAR**

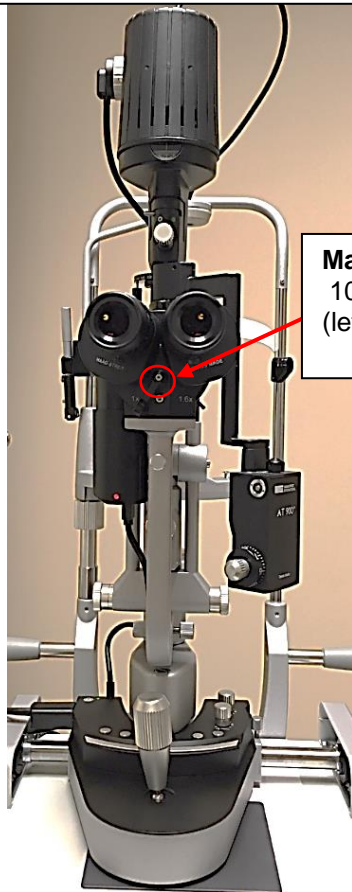

**Magnification:**  
10X or 16X  
(left) (right)

For optimal views,  
candidates must use  
focus rod to focus the  
oculars

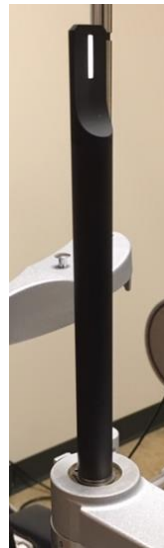

### Bright White:


Top pin fixed all the  
way to the left



### Power:


blinking green = **standby**  
solid green = **on**

**Non-functioning buttons**

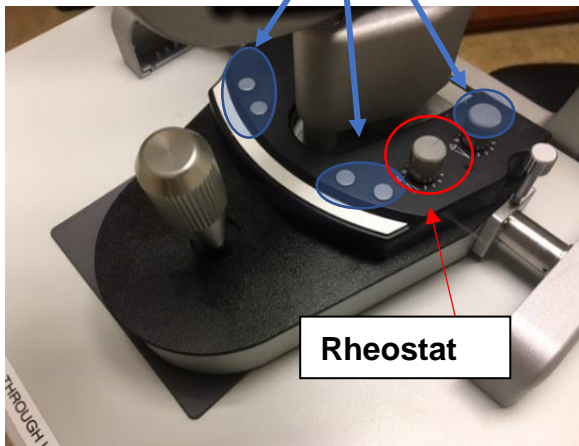

**Rheostat**

## OR STATIONS SLIT LAMP OVERVIEW

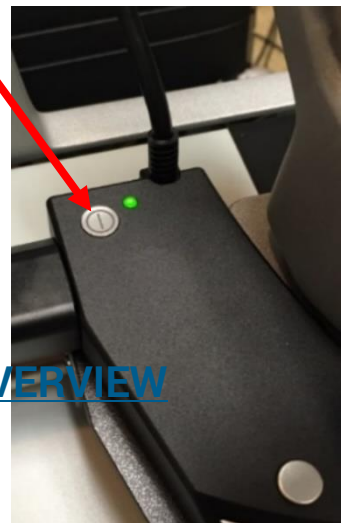



### **Bright White**

*Largest Aperture - No Filters*

Top pin pushed to the **left** all the way



### **Blue Filter**


Circular dial turned all the way to the **left** so "B"  
is displayed

### **Yellow Filter**

If the lever is at , the filter is IN

If the lever is at , the filter is OUT

(the filter cannot be in between these positions)

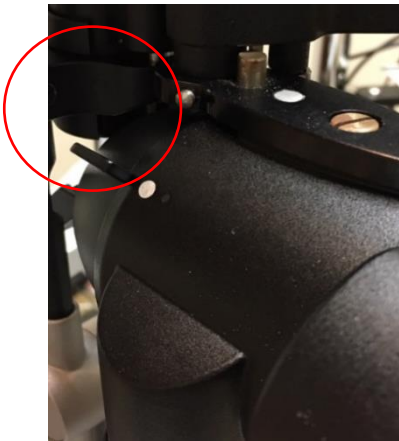



## ANTERIOR STATION EQUIPMENT & SUPPLIES






## ANTERIOR STATION EQUIPMENT & SUPPLIES

**Note: The brand of the item may change and be different than what is displayed**

Example of anterior segment station supplies set up on the counter ready for use

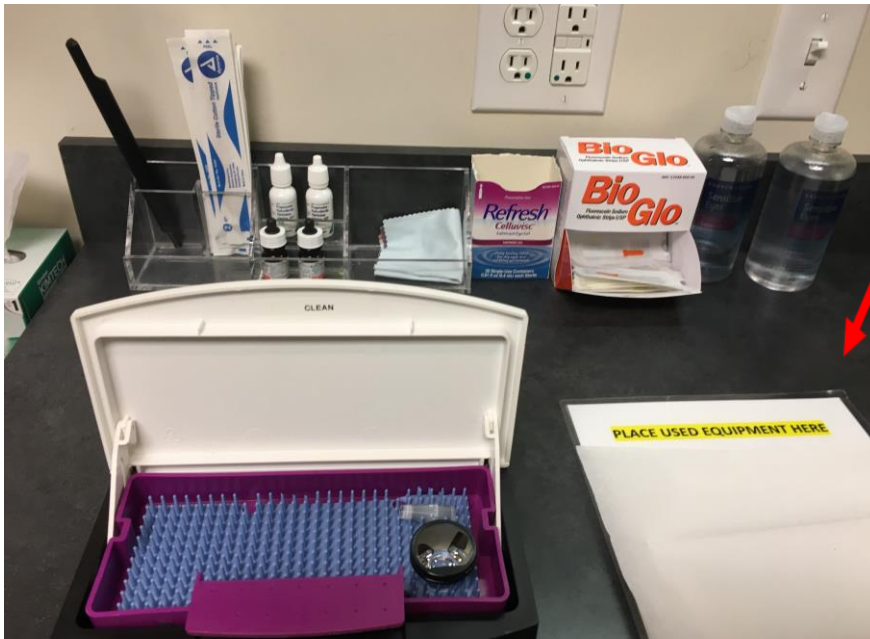

### Tray

- Tray will contain ready to use, already-disinfected, gonio lens and tonometer
- Please keep tray closed until the start of the exam
- After using the instruments please place on the clean paper towel



## ANTERIOR STATION EQUIPMENT & SUPPLIES

**Note: The brand of the item may change and be different than what is displayed**

### **Gonio Lens**

Volk 4-Mirror Glass Gonio  
Fundus Lens with Flange



### **Kim Wipes**



### **Elbow Rest**




### **Alcohol Wipes**



### **Repeat Policy**

Posted on the back of the **skills station** exam room doors, next to the computer workstation

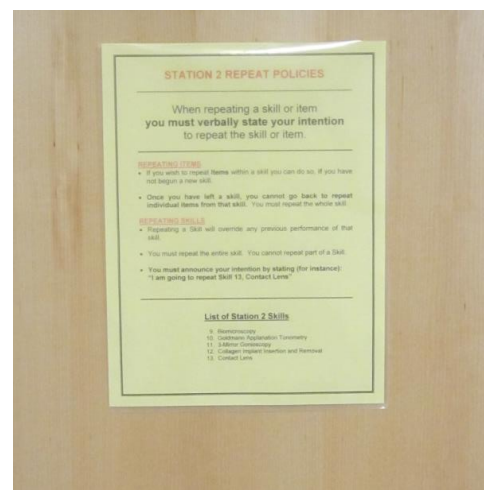

## POSTERIOR STATION EQUIPMENT & SUPPLIES





## POSTERIOR STATION EQUIPMENT & SUPPLIES

### **"Modified" Keeler BIO Vantage Plus Slimline Wireless Plus Indirect System**



- The BIO headset will be hanging on the stand
- Without an appropriate target and working distance (18-25 inches) the target will appear decentered

#### **FIXED**

- aperture
- filter
- brightness
- circle beam height
- +2.00 oculars

#### **ADJUSTABLE**

- fit behind head
- fit on top of head
- PD
- Pantoscopic tilt: angle of ocular to eyes


Illumination  
will be set by  
the SP

Slide to adjust  
PD between  
48-76

Turn to tighten or  
loosen fit

Knob loosens to adjust  
ocular distance to eyes

Vertical fit

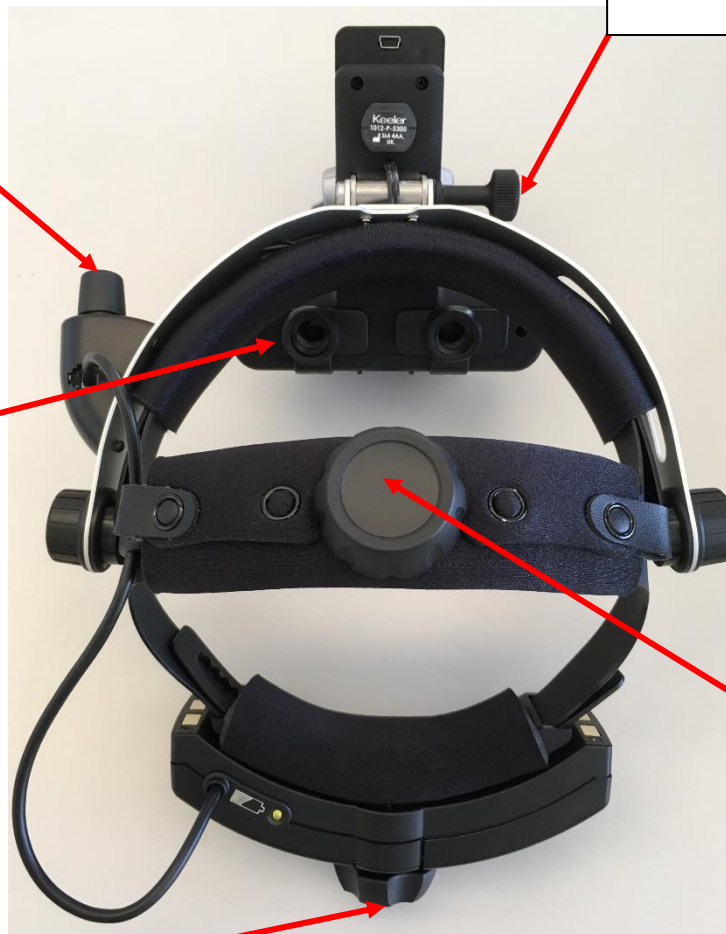

# POSTERIOR SEGMENT STATION EQUIPMENT AND SUPPLIES

## Lens Overview

The following lenses are available for candidate use during posterior evaluation

- Volk 20D Lens
- Volk VDGTL1 Digital High Mag Imaging Lens
- Volk Digital Wide Field Lens
- Volk 78D & 90D Lenses



## Miscellaneous Supplies

Alcohol Wipes



Kim Wipes



Tissues



Elbow Rest



## Repeat Policy

Posted on the back of the **skills station** exam room doors, next to the computer workstation

